
# Operation Microscope OMS-710





**CONNECTING VISIONS** 

# Excellence in image and operability Operation mircroscope OMS-710

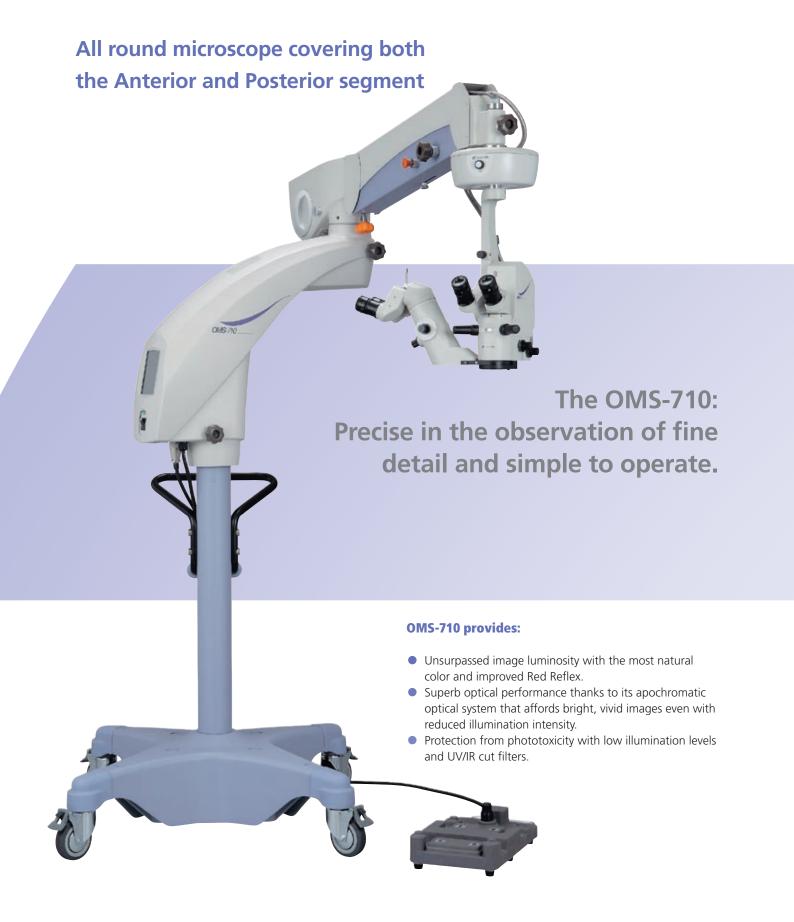


## **Improved Red Reflex and Shadow Contrast**

The newly designed optics and illumination system offer unsurpassed Red Reflex quality.

With combined 2° and 4° illumination, contrast is more natural and the new color temperature in the illumination system, allows a softer and more detailed view.



Red reflex with Color temperature filter

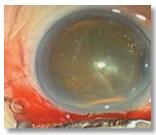

View without Color temperature filter



# TOPCON

The OMS-710 incorporates Illumination Shutter on 4° illumination. Preferable shutter position can be set to enhance Red Reflex.

## **Features**

## **Main Microscope Binoculars**

Two types of binoculars are available for the OMS-710:

- 1. A 45° fixed set of binoculars
- 2. A variable inclination set of binoculars; adjustable from 45° to 90°, providing a more comfortable observation image

## **Compact Body**

The OMS-710's compact base and longer arms allow easy mobility and storage in limited space.

## **XY Translator**

The precise and smooth movement of the X/Y translator provides the surgeon with control and ease of use. The centering button lights up when zoom, focusing mechanism and X/Y movement are set to the initial position.




## **Touch Panel LCD Monitor**

A sophisticated touchscreen LCD monitor allows for a quick and easy change of settings and programs.



# **Optional Accessories**

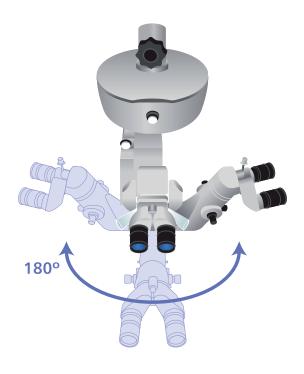

## **Rotating Assistant Microscope**

The Assistant Microscope easily swings from one side to the other, pivoting comfort, flexibility and convenience.

## **Compact TV Relay Lens**

The compact TV Relay Lens allows to connect 1/2 inch and 1/3 inch CCD cameras either with a C-mount or Sony-2 mount blocking or interfering with the Assistant's microscope movement. When a CCd camera is mounted on the OMS-710, there is no interference with any functionality.


## **Coaxial Binocular Observation Tube**


The coaxial binocular tube provides the same view as the observation tube of the surgeon. This is useful for education and teaching purposes. Easy to position and simple to adjust, the observation height is identical to the surgeon's.

## **Diaphragm Changer**


The flip-in diaphragm provides a wider focusing range. Two steps of diaphragm are available.

# **Dimensions**

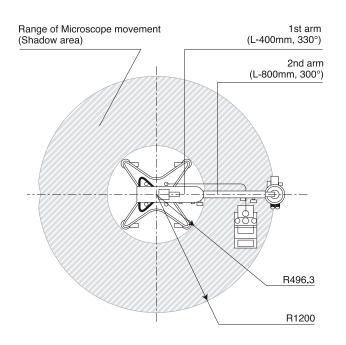

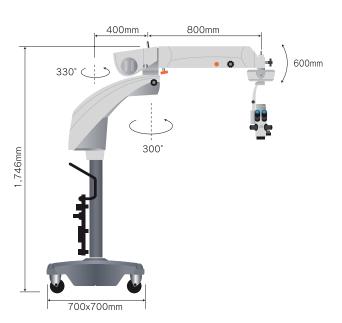

# **System Chart**

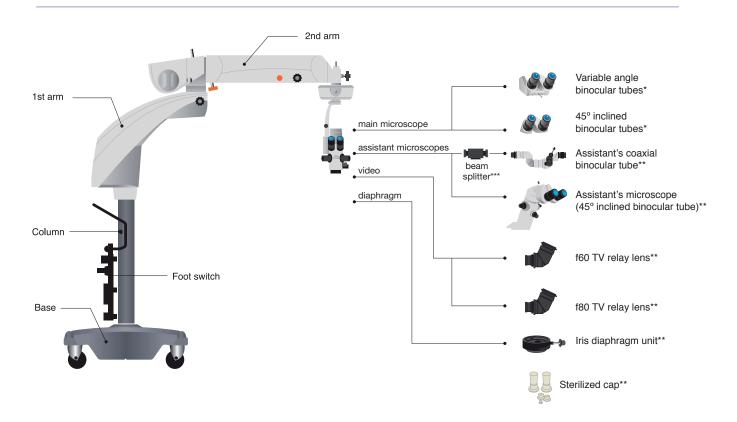

- \*To be selected at purchase

  \*\*Optional accessories

  \*\*\*Included with variable angle binocular type as standard

## **Specifications**

Eyepiece	12.5x
Objective Lens	f=200 mm
Magnification Selection	Continuous zoom
Magnification Indication	4.2/5/6/7/8/9/10/11/13/15/17/19/21
Total Magnification	4.2-21x
1st Arm Length / Rotation	400mm / 330°
2nd Arm Length / Rotation	800mm / 300°
Arm Load	9.5-16.0 kg
Accessory Capacity	5.5 kg
Power Voltage	AC100-240V 50/60Hz
Filter	Built-in: Heat-absorbing and UV cut filters
	Flip-in/out: Yellow and Color temperature conversion filters
Size	Base: 700 x 700 mm
	1,746 mm
Weight	45° Inclined Binocular Type: 184Kg
	Angle Binocular Type: 185Kg





Subject to change in design and/or specifications without advanced notice.

In order to obtain the best results with this instrument, please be sure to review all user instructions prior to operation.

## Topcon Europe Medical B.V.

Essebaan 11; 2908 LJ Capelle a/d IJssel; P.O. Box 145; 2900 AC Capelle a/d IJssel; The Netherlands Phone: +31-(0)10-4585077; Fax: +31-(0)10-4585045 E-mail: medical@topcon.nl; www.topcon.eu

### **Topcon Danmark**

Praestemarksvej 25; 4000 Roskilde, Danmark Phone: +45-46-327500; Fax: +45-46-327555 E-mail: topcon@topcondanmark.dk www.topcondanmark.dk

## Topcon Scandinavia A.B.

Neongatan 2; P.O. Box 25; 43151 Mölndal, Sweden Phone: +46-(0)31-7109200; Fax: +46-(0)31-7109249 E-mail: medical@topcon.se; www.topcon.se

### Topcon España S.A.

HEAD OFFICE; Frederic Mompou, 5; 08960 Sant Just Desvern; Barcelona, Spain Phone: +34-93-4734057; Fax: +34-93-4733932 E-mail: medica@topcon.es; www.topcon.es

## Topcon Italy

Viale dell' Industria 60; 20037 PADERNO DUGNANO; (MI) ITALY Phone: +39-02-9186671; Fax: +39-02-91081091 E-mail: topconitaly@tiscali.it; www.topcon.it

## Topcon Deutschland GmbH

Giesserallee 31; D-47877 Willich, Germany Phone: +49-(0)2154-8850; Fax: +49-(0)2154-885177 E-mail: info@topcon.de; www.topcon.de

#### Topcon (Great Britain) Ltd.

Topcon House; Kennet Side; Bone Lane; Newbury Berkshire RG14 5PX; United Kingdom Phone: +44-(0)1635-551120; Fax: +44-(0)1635-551170 E-mail: info@topcon.co.uk; www.topcon.co.uk

## Topcon Ireland

Unit 276, Blanchardstown; Corporate Park 2 Ballycoolin; Dublin 15 Phone: +353-18975900; Fax: '+353-18293915 Email: medical@topcon.ie; www.topcon.ie

## Topcon S.A.R.L.

HEAD OFFICE; 89, rue de Paris; 92585 Clichy, France Phone: +33-(0)1-41069494; Fax: +33-(0)1-47390251 E-mail: topcon@topcon.fr; www.topcon.fr

# Topcon Portugal

Rua da Forte, 6-6A, L-0.22; 2790-072 Carnaxide; Portugal Phone: +351-210-994626 Fax: +351-210-938786 www.topcon.pt

### Topcon Polska Sp. z o.o.

ul. Warszawska 23; 42-470 Siewierz; Poland Phone: +48-(0)32-670-50-45; Fax: +48-(0)32-671-34-05 www.topcon-polska.pl

