
SOLOS Series

Automatic Lens Analyser

Drive Premium Lens Sales with Advanced Automated Analysis





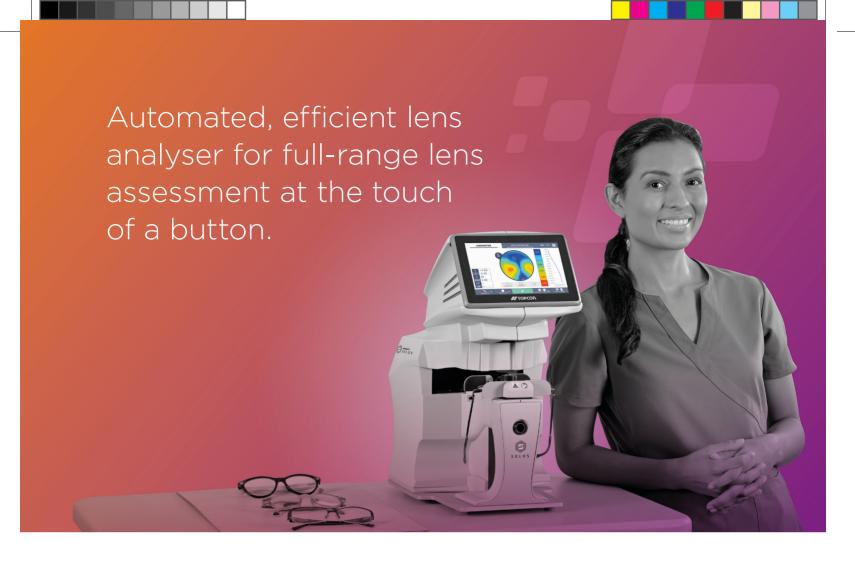

VERSATILITY


- Seamlessly detect, measure and mark single vision, progressive, and multifocal* lenses.
- Detailed assessment of both glazed and uncut lenses, catering to a range of clinical and optical needs.

ADVANCED LENS INSIGHTS


- Comprehensive lens analysis capabilities, including mapping optical features (e.g., distortion areas) and generating comprehensive aberration maps for progressive lenses.
- Precise evaluation of light absorption and colour appearance with a full-range spectrophotometer.**
- Enhanced understanding of lens performance and advantages of premium lens designs.

EFFICIENCY

- Intuitive, easy-to-use interface delivering efficient operation and rapid onboarding.
- One-touch functionality enables full delegation.
- Seamless data sharing to EMRs, the CV-5000PRO digital refraction system, or Chronos binocular refraction system.
- Detailed PDF reports to support premium lens discussions at dispensing workstations.

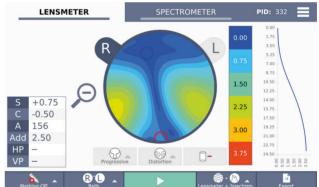

Detailed Lens Mapping

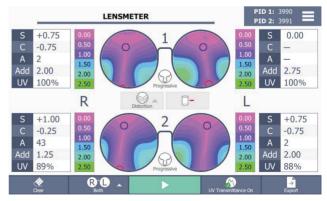

Comparison Environment (screenshot from SOLOS Lite)

^{*} SOLOS and SOLOS Lite are not compatible with trifocal lenses

^{**} On the SOLOS Lite, only UVA transmittance at a single wavelength is measured (no full-range spectrometer and chromaticity data available)

OVERVIEW



Automated, One-Touch Operation



Lens Mapping with **Distortion Map**



UVA, Blue Light, and Visible Light **Transmittance** Measurements**



Automatic Lens Marking



Automatic Lens Type Detection




WirelessData Transfer



Extended **Measurement** Range (Up to +/- 20D)




Detailed **PDF** Reports

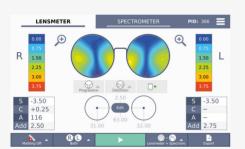


Lens Refractive Power Analysis

Featuring

- Automatic detection of optical centres
- Distortion maps highlighting progressive and free-form lens features
- Option to adjust the pupillary distance horizontally to view prisms
- Toggle between lensmeter and spectrometer data**

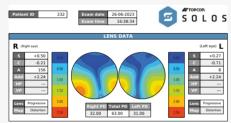



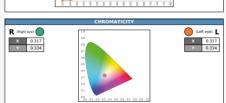

Spectacle Lens PDF Report

Features and benefits

- Lens Data
 - Lens type
 - Distortion map
 - Monocular and binocular lens centration distances
- Transmittance**
 - UVA, Blue and Visible Light Transmittance
 - Illustrate adequate protection against harmful UV rays
 - Showcase lens performance for various wavelengths for both prescription lenses and sunglasses
 - Support discussion of blue light blocking lenses
- $\bullet \ \, \text{Chromaticity}{}^{**}$
 - Lens colour specified using internationally recognised CIE chromaticity coordinates
 - Verify colour consistency for tinted lenses

INTRODUCING SOLOS LITE


Advanced lens analysis with exact UVA transmittance at a single wavelength. Experience the trusted reliability of SOLOS with our entry-level model.



FEATURES	SOLOS	SOLOS Lite
Fully Automatic Lens Analysis	Yes	Yes
Automatic Lens Type Detection	Yes	Yes
Automatic Pupil Distance Measurement	Yes	Yes
Automatic Lens Marking	Yes	Yes
Comparison maps	Yes	Yes
UVA spectrometer	Yes	Yes
Full range spectrometer	Yes	No

GENERAL DATA			
Dimensions	W245 mm x H450 mm x D354 mm		
Weight	Net Weight: 8.0 kg, Cables and user manual: 2.0 kg, Packaging: 3.2 kg		
Printer	Internal (thermal)		
Screen	7" Touch screen LCD/16M colours		
Light source	Green e-line source		
Working conditions	10°C to 40°C		
Power supply	AC 100 - 240 V — 50/60 Hz		
Classification	Class I Medical Device - EU Regulation 2017/745 and US FDA 21 CFR		
Standards	ISO 8598		
Data output	LAN, Wi-Fi		
MEASUREMENT RANGE			
Sphere power	-20D — +20D (step 0.01, 0.0625, 0.125, 0.25D)		
Cylinder power	-10D — +10D (step 0.01, 0.0625, 0.125, 0.25D)		
Cylinder axis	0 - 180° (step 1°)		
Addition power	-4D — +4D (step 0.01, 0.0625, 0.125, 0.25D)		
Prism power	OD - +20D (step 0.01, 0.0625, 0.125, 0.25D)		
PD measurement	Mono / Bino		
Cylinder notation	- / +		
OTHER			
Spectrometer	SOLOS 315nm - 800nm	Solos Lite 365nm (single wavelength and no chromaticity)	
Automatic marking	Optical centre and axis (framed spectacles or uncut lenses)		
MEASURED LENS SIZE			
Single lens	Diameter: 22mm —80mm	Thickness: up to 22mm	
Spectacles	Height: 22mm —80mm	Wrapping angle: up to 18°	

CONNECTIVITY

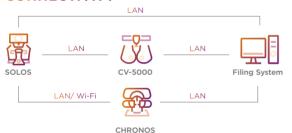





TOPCON HEALTHCARE UNIVERSITY

Eye Health Education Begins Here: learning.topcon.com

Note: Subject to change in design and/or specifications without advanced notice

TOPCON CORPORATION75-1 Hasunuma-cho, Itabashi-ku, Tokyo 174-8580, JAPAN.
Phone: +81-(0)3-3558-2522/2502 Fax: +81-(0)3-3965-6898 www.topconhealthcare.jp

TOPCON HEALTHCARE SOLUTIONS EMEA OY

Saaristonkatu 9, 90100 Oulu, FINLAND Phone: +358-20-734-8190 E-mail: thsemea.sales@topcon.com www.topconhealthcare.eu

TOPCON EUROPE MEDICAL B.V.

Essebaan 11, 2908 LJ Capelle a/d IJssel THE NETHERLANDS Phone: +31 -(0)10-4585077 Fax: +31 -(0)10-4585045 E-mail: medical@topcon.com www.topconhealthcare.eu

TOPCON DANMARK

Praestemarksvej 25, 4000 Roskilde DANMARK Phone: +45-46-327500 Fax: +45-46-327555 E-mail: info.todk@topcon.com www.topconhealthcare.eu/dk_DK

TOPCON SCANDINAVIA

Neongatan 2, P.O.Box 25, 43151 Mölndal SWEDEN Phone: +46-(0)31-7109200 Fax: +46-(0)31-7109249 E-mail: info.hcs.se@topcon.com www.topconhealthcare.eu/sv SE

TOPCON ESPAÑA

Oficina Principal en España Frederic Mompou, 4, 08960 Sant Just Desvern, Barcelona, SPAIN Phone: +34-93-4734057 Fax: +34-93-4733932 E-mail: medica@topcon.com www.topconhealthcare.eu/es_ES

Viale dell' Industria 60, 20037 Paderno Dugnano, (MI) ITALY Phone: +39-02-9186671 E-mail: info_tit@topcon.com www.topconhealthcare.eu/it_IT

VISIA IMAGING S.R.L.
Via Martiri della Libertà 95/e, 52027 San Giovanni Valdarno (AR), ITALY

TOPCON FRANCE MEDICAL

1 rue des Vergers, Parc Swen, Bâtiment 2, 69760 Limonest, FRANCE Phone: +33-(0)4-37 58 19 40 Fax: +33-(0)4-72 23 86 60 E-mail: topconfrance@topcon.com www.topconhealthcare.eu/fr FR

TOPCON DEUTSCHLAND MEDICAL

Hanns-Martin-Schleyer-Strasse 41, D-47877 Willich, GERMANY Phone: (+49)2154-885-0 Fax: (+49)2154-885-177 E-mail: info@topcon-medical.de www.topconhealthcare.eu/de_DE

TOPCON POLSKA SP. Z. O. O.

ul. Warszawska 23, 42-470 Siewierz POLAND Phone: +48-(0)32-670-50-45 Fax: +48-(0)32-671-34-05 E-mail: info.tpl@topcon.com www.topconhealthcare.eu/pl_PL

TOPCON GREAT BRITAIN MEDICAL

Topcon House, Kennet Side, Bone Lane, Newbury, Berkshire RG14 5PX UNITED KINGDOM Phone: +44-(0)1635-551120 Fax: +44-(0)1635-551170 E-mail: medical.tgbm@topcon.com www.topconhealthcare.eu/en_UK

TOPCON IRELAND MEDICAL

Unit 292, Block G, Blanchardstown, Corporate Park 2 Ballycoolin Dublin 15, D15 DX58, IRELAND Phone: +353-12233280 E-mail: medical.ie@topcon.com www.topconhealthcare.eu/en_IE